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USER GUIDE

CERAGEM MEDISPA PRO Skin Care Device
USER MANUAL

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CERAGEM MEDISPA PRO USER MANUAL

- This device is a therapeutic massager for relaxation for the face and also provides skincare.
- In order to use the product safely and correctly, be sure to read this user manual before using it.
- This manual includes the product warranty.
- In the future, we may improve upon the design and specifications of this product without notice.
- "CAUTION" and "Risk of electric shock, do not remove cover (or back). Refer servicing to qualified service personnel".

	Usage and Storage Caution	04
01	Precautions for Use and Storage	06
02	Name of Each Part and Description of Functions ..	10
03	How to Use the Product	18
04	Troubleshooting	28
05	FAQ	29
06	Product Specification	30
07	Product Warranty	31

Usage and Storage Caution

Do not apply any skincare product onto the target skin area before using this device!

Please read through the Directions for Use carefully and keep it in a safe place for future reference. You may contact us online if you have any questions during the use of this device.

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following :

READ ALL INSTRUCTIONS BEFORE USING, SAVE THESE INSTRUCTIONS

Danger

In cases of violating the instructions, there is an emergent possibility of serious damage or even death.

Warning

In cases of violating the instructions, there is some possibility of serious injury or death.

Attention

In cases of violating the instructions, there is possibility of minor injury or property damage.

Danger

To reduce the risk of electrocution :

- (1) Always unplug this product immediately after using it.
- (2) Do not use while bathing.
- (3) Do not try to insert any object into this device.
Refrain from leaving this device untended anywhere.
- (4) Do not use this device around the eyes.
- (5) Due to a lack of data regarding the use of this device on children aged below 15 years, it is advised to keep this device out of touch by children aged below 15 years. Similarly, please use this device at a place far from the reach of these children.
- (6) Do not use this product near any water-filled bathtub/basin, washbasin or other similar utensils.
- (7) Do not place or store product where it can fall or be pulled into a tub or sink.
- (8) Do not place in or drop into water or other liquid.
- (9) Do not reach for a product that has fallen into water. Unplug immediately.

Warning

To reduce the risk burns, electrocution, fire, or injury to persons.

- (1) A product should never be left unattached when plugged in.
- (2) Close supervision is necessary when this product is used by, on, or near children or invalids.
- (3) Use this product only for its intended use as described in this manual.
Do not use attachments not recommended by the manufacturer.
- (4) Never operate this product if it has a damaged cord or plug, if it is not working properly, if this has been dropped or damaged into water.
Return the product to a service center for examination and repair.

- (5) Keep the cord away from the heated surface.
- (6) Never block the air opening of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- (7) Never use it while sleeping or drowsy.
- (8) Never drop or insert any object into any opening or hose.
- (9) Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- (10) Connect this product to a properly grounded outlet only. See Grounding Instructions.
- (11) Unplug this product before filling. Fill (reservoir) with water only unless otherwise specified by manufacturer. Do not overfill (or specify filling instructions).
- (12) Please choose an appropriate level with your skin conditions. When using this product for the first time, please start it up from Level 1.
- (13) After being started up, this product shall not be left untended. Be sure to turn this product off in time after using it.
- (14) Before moving this product from one place to another where the temperature is higher or lower by up to 68 degrees Fahrenheit (e.g. from a low-temperature place to a high-temperature place, or vice versa), turn it off and leave it unused for 2 hours or so before starting it up again for re-use.
- (15) Do not use this product near any flammable or explosive substance (inclusive of any sprayer or petrol station, etc.) to prevent bodily injuries caused by the device induced fire accident or electric shock.
- (16) Please check the room temperature before using this product which shall be used in a place with an ambient temperature of 32-104 degrees Fahrenheit.
- (17) Do not expose the target skin area to ultraviolet radiation directly within 48 hours after hair removal.

- (18) Do not expose this product to any intensive sunlight or a functioning ultraviolet lamp for long hours which will do severe harm to it.
- (19) The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.
- (20) Children being supervised not to play with the appliance.
- (21) For a part of class III construction supplied from a detachable power supply unit, the instructions state that the appliance is only to be used with the unit provided.

Attention

Cleaning and Maintenance

In order to ensure this device to deliver its best performance and ensure the longest possible service life, be sure to clean up this device before and after using it every time so as to prevent dust from accumulating on its device for long hours. This device would be unable to deliver its promised performance should you fail to clean it up.

- (1) Examine this device, its head and component parts of its body in particular and ensure that all of them remain clean and tidy.
- (2) Do not rinse out the device or any of its fittings with running tap water, in a shower, or inside a dishwasher.
- (3) Do not use a steel wire sponge, any abrasive detergent (e.g. toothpaste, etc.) or corrosive liquid (e.g. gasoline, engine oil or acetone) to clean up this device.
- (4) Switch off and unplug this device and let it cool down. When this device is being used, its light outlet will become warm. Be sure to let this device cool down before cleaning it up.
- (5) Where necessary, use a soft piece of cleaning cloth to clean up the exterior of this device.
- (6) Dab a few drops of water onto the cleaning cloth and then use it to clean up the following component parts: the light outlet and the cooling panel.

01. Precautions for Use and Storage

Before use

Any of the following people should take caution when using the product.

- (1) For people with skin disease, who have health conditions such as allergies or sensitive skin, or who are receiving skin treatments due to skin-related procedures, etc., it is recommended that they consult a specialist before using the product.
- (2) For people who have undergone plastic surgery or skin-related procedures or treatments, it is recommended that they use this product after the skin has recovered (approximately 2 months).
- (3) It is not recommended that people who are pregnant, lactating or who have weak immune systems use this product.
- (4) It is not recommended that people who wear or have a cardiac pacemaker implanted use this product.
- (5) Do not allow infants and children under the age of 12 to use this product.
- (6) Regardless whether or not you use the product, your skin condition may change temporarily due to seasonal changes, external environmental factors, alcohol consumption, and physical factors such as stress. If your skin condition has changed, it is recommended that you stop using this product until your skin has returned to its normal condition or after consulting a specialist.
- (7) Stop using this product if your face turns red, or repeatedly has itching or purulent inflammation while using the product.
- (8) Do not use this product immediately after shaving, hair removal, or on wounds or other areas with abnormal skin reactions.
- (9) It is recommended that you observe the recommended time of use.
- (10) Before performing skin care, remove any jewelry from the areas this product will come into contact with.

- (11) Do not use this product while in motion, such as driving, walking, etc.
- (12) This product is a skin care device. Do not use it for purposes other than its original purpose.

When using

Any of the following people should take caution when using the product.

- (1) If you fall under any of the following categories, do not use this product.
 - Persons with metal allergies
 - Persons with nociceptive or perceptual disturbance
 - Persons with severe heart disease
 - Persons with a high fever
 - Persons with severe abnormal blood pressure
 - Persons who have used steroid hormones for a long period of time
 - Persons with epilepsy
 - Persons with sensitive skin or allergic skin
 - Persons who have had dental implants or who are using metal prostheses
 - Children under the age of 12
 - Pregnant women
 - Persons with telangiectasis
 - Persons with infectious diseases
 - Persons with other skin diseases

- (2) Do not use this product on any of the following areas.
 - Areas with bleeding
 - Areas with severe acne or skin problems
 - Ears, eyelids, throat and mucous areas (mouth, inside the nose, genitals, ocular mucosa, etc.)
 - Areas with a wound
 - Areas around bone fractures
 - Areas with various severe skin problems
 - Areas that have an electronic device implanted

How to use and handle

- (1) It is recommended that this product be used once a day, and that the instructions are followed when used.
- (2) Be careful not to use this product in one area for a long period of time.
- (3) Do not use this product immediately after any procedure or surgery, and consult a specialist before use.
- (4) Avoid using this product on areas around your eyes, Adam's apple and philtrum.
- (5) Gaps on the surface that comes into contact with the skin are waterproof, but other gaps on the product are not waterproof, so be careful not to allow cosmetic products or water to enter the product.
- (6) Pressing the buttons too hard may cause a failure.
- (7) Do not attempt to disassemble, repair or modify this product. Doing so may cause fire, electric shock or injury due to malfunction.
- (8) Take caution while handling this product as dropping it or applying impact to it may cause failure or damage.
- (9) After using the product, if there is a foreign substance on the area that comes into contact with the skin (head), use a soft cloth, gauze, wet wipes, cotton pads, cotton swabs, etc., to wipe the product. Do not wipe the product surface too forcefully and do not use strong chemical products such as alcohol, benzene or thinner, etc., as doing so may alter the product surface.
- (10) Store the head unit in clean and dry conditions without any foreign objects.
- (11) Be careful not to have the product be exposed to direct sunlight in places such as inside a closed car.

Use and Storage Conditions

Storage and Transporting Conditions	Operating Temperature	32°F - 95°F
	Humidity	20% - 80%RH
	Relative Pressure	700hPa - 1060hPa
Operating Conditions	Storage Temperature	-4°F - 122°F
	Humidity	-95% RH
	Relative Pressure	700hPa - 1060hPa

Precautions For Storage

Do not store the product in any of the following places. Doing so may cause fire, electric shock or failure.


- Places with direct sunlight or where heat is generated such as on an electric heater, etc.
- Places where the temperature is too high (over 104°F) or too low (below 32°F)
- Places with excessive moisture, dust or soot, etc.
- Places near flammable substances or volatile liquids
- Places where the product can be immersed in water or other liquids such as a drain or sink, etc.
- Dishwasher, oven, washing machine, etc.
- Places with severe vibrations

Disposing of the Product


When disposing of the product, you must recycle it in a manner that complies with the waste electronics (WEEE) guidelines. When disposing of the product, dispose of the main body and the cable separately in accordance with the regulations.

02. Name of Each Part and Description of Functions


1) Main Body

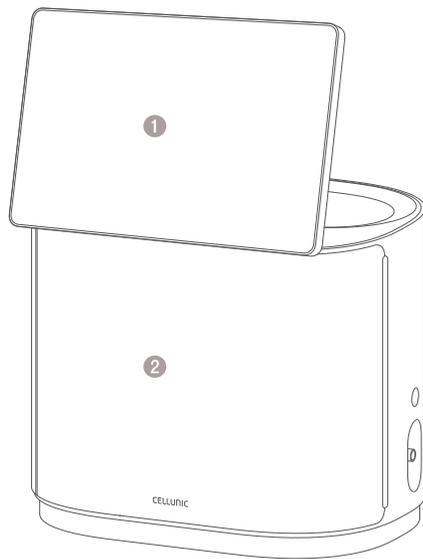Front of the Main Body

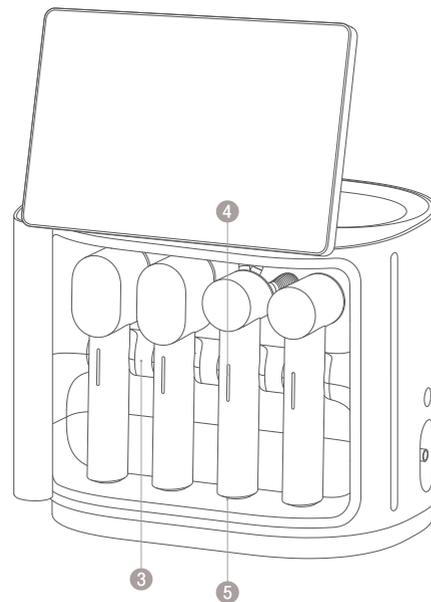Inside the front of the main body

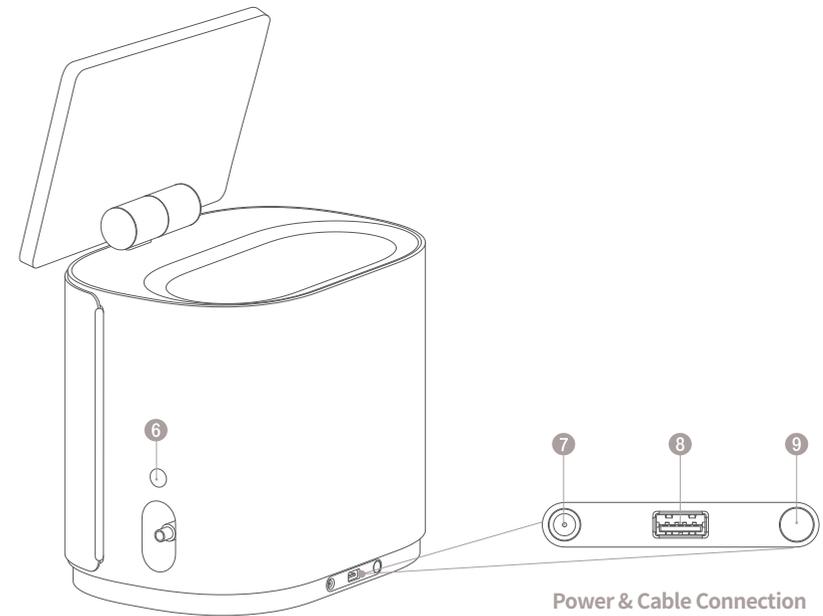
Side/back of the main body

Name	Functions
① LCD Screen	Select among Personal mode and Expert mode, set and run Care
② Main Body Cover	Cover for protecting the handpiece holder, used to store and prevent contamination of the handpiece
③ Handpiece Holder	Holder for storing the handpiece
④ Handpiece LED	A white LED turns on when the handpiece is connected correctly A red LED turns on when the handpiece is connected incorrectly
⑤ Holder LED	When the handpiece is used, an LED turns on in accordance with the order of care


Name	Functions
⑥ Handpiece Cable & Returner	The cable that is connected to the handpiece and the returner button for rewinding the cable
⑦ Adapter Connector	The connector that is connected to the adapter that supplies power
⑧ USB Port	A port for sending usage history and updating the product by connecting a USB device
⑨ LCD Power Button	You can enter or cancel screen sleep by pressing and holding it down for more than 2 seconds.

2) Handpiece and Components

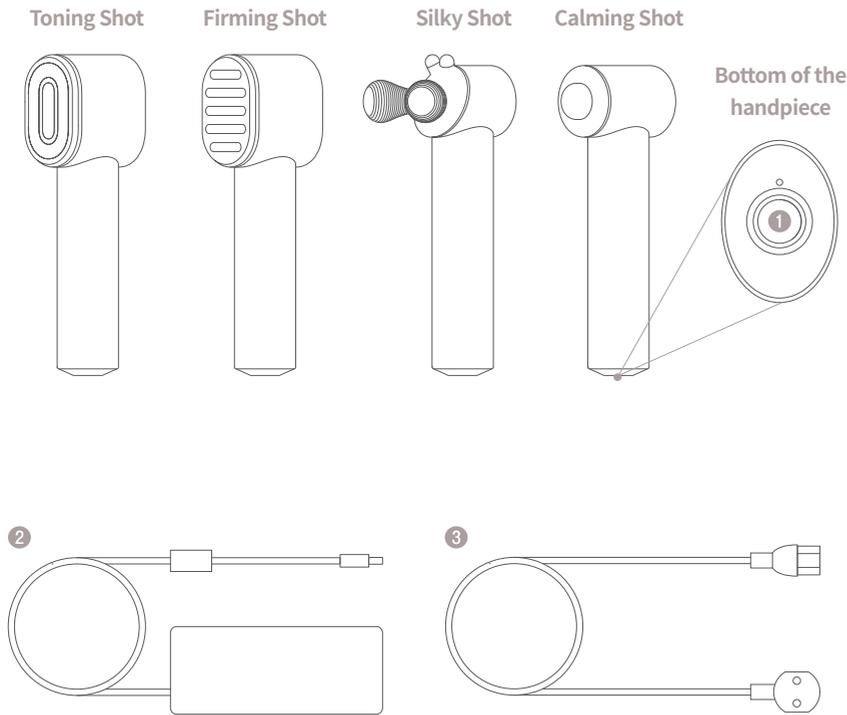
Name	Functions
1 Cable Connector	A connector that connects each handpiece to the cable
2 Adapter	An adapter for supplying power to the main body
3 Power Cord	A cord that connects the power to the adapter

3) Description of Functions General Mode

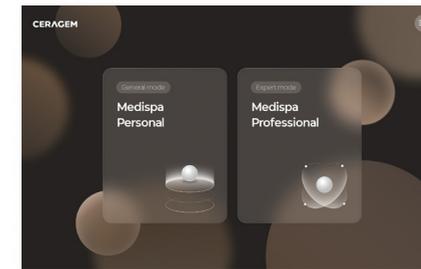
(1) Home Screen

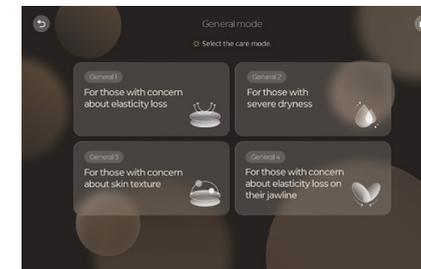(2) General Mode Start Screen

4 modes selection screen

- For those concerned with a loss of elasticity
- For those with severe dryness
- For those concerned about skin texture
- For those concerned with a loss of elasticity of the jawline

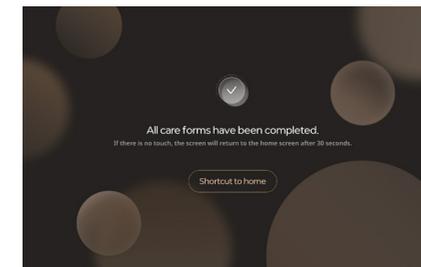
(3) Ending Care

If you do not touch the screen after care in general mode has ended, the screen will return to the home screen after 30 seconds.

3) Description of Functions Expert Mode


(1) Home Screen

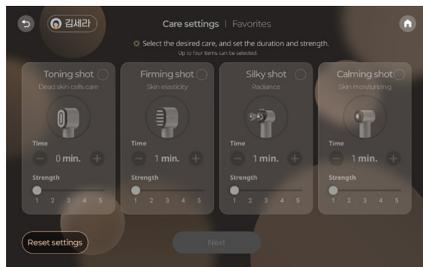(2) Expert Mode Start Screen

- The desired time (up to 10 minutes) can be set
- The desired strength (up to 5) can be set

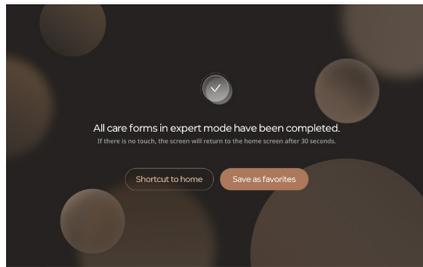(3) Ending Care

If you do not touch the screen after care in expert mode has ended, the screen will return to the home screen after 30 seconds.

3) Description of Functions LCD screen


(1) Sleep Mode

When the Power button at the bottom of the main body is pressed and held down for 2 to 3 seconds, the mode will switch to the sleep mode or the sleep mode will be canceled.

(The LCD does not operate and the main body switches to stand-by status)

- ◀ Sleep mode

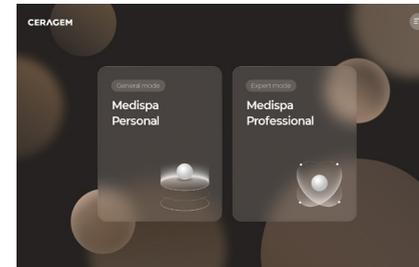- ◀ When canceling sleep mode


(2) Screen Save Mode

If the device does not operate for 5 minutes, the mode will switch to Screen Save mode.

When you knock on the screen once, Screen Save mode will be canceled.

(The LCD & the main body switch to stand-by status)

- ◀ Screen Save mode

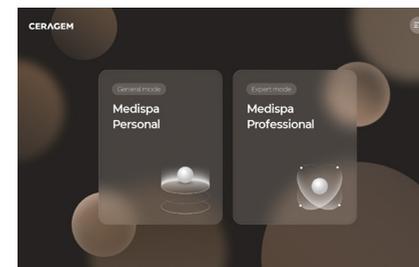- ◀ When Screen Save mode is canceled


4) Description of the Main Screen UI

Care In Progress Screen

1	Home button	Ends the care in progress and returns to the Home screen.
2	Video Mode Switch button	During operation, you can check the correct instructions for using the product by pressing the 'Switch to Video Mode' button.
3	Care Timer	You can check the remaining operating time of the current handpiece on the timer.
4	Current Care Step	You can check which step the currently operating handpiece is in among all of the care steps.
5	Handpiece Status	You can check the operating status of the handpiece. ON is displayed when the handpiece is turned on, and OFF is displayed when it is turned off.
6	Strength Adjustment Bar	After pausing the care, you can change the strength of the handpiece. The strength cannot be changed while the handpiece is operating.
7	Pause/Continue button	After starting operation, when you press Pause mid-operation to pause the handpiece and press the 'Continue' button displayed where the 'Pause' button is, the relevant care will resume for the remaining time. When pausing, the handpiece LED will flash.
8	Skip button	When you press the 'Skip' button, the screen will move to the Connect to the Next Handpiece screen.

Expert Mode Care Settings Screen

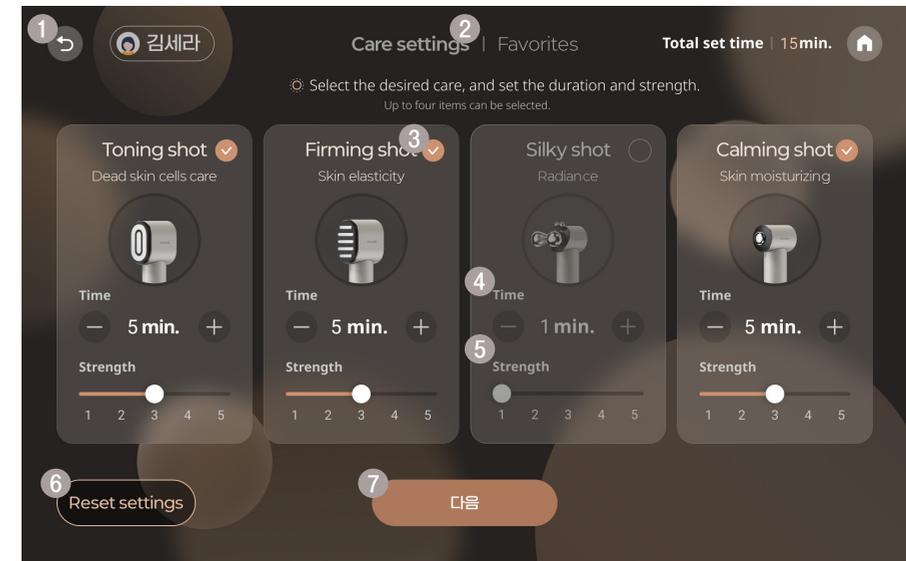1	Previous Button	The screen returns to the Login screen, which is the previous screen.
2	Care Settings/Favorites Button	From the 'Care Settings' tab, you can set and carry out the desired details of the care, The 'Favorites' tab cannot be used while using the product in Guest mode, and up to 3 details of the care you set can be saved.
3	Handpiece Selection Button	You can select the desired handpiece by pressing the circle. When you press it again after the handpiece has been selected and marked with a check mark, the handpiece will be deselected.
4	Time Adjustment +/- Button	You can select how many minutes you will operate the handpiece for. (Up to 10 minutes)
5	Strength Adjustment Bar	You can set the operating strength of the handpiece. The operating strength can also be changed after care has started.
6	Initialize Settings Button	When you initialize settings, all the selected handpieces will be deselected and time and strength settings will be initialized.
7	Next button	The screen will move to the screen where you can check the details of the care you set. ※ The design of the screen is subject to change without prior notice for performance improvements.


03. How to Use the Product

1) Turning the Product On

(1) Connecting the Adapter

When the adapter is connected to the main body, the product will be turned on.

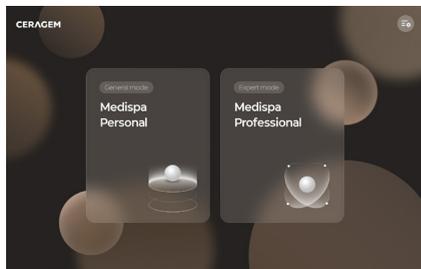(2) Selecting the Mode


After the power is turned on, the initial screen will appear. From the initial screen, you can select and execute the desired function among 'General mode' and 'Expert mode'.

2) How to Carry out Care in General Mode

(1) Selecting General Mode

From the Home screen, select General mode.

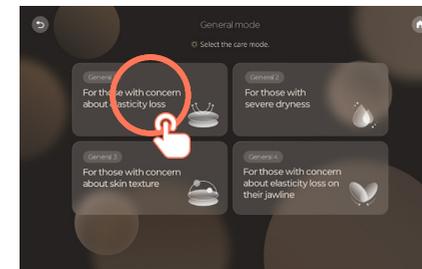
(2) Selecting the Desired Care Mode

In General mode, there are a total of 4 Care modes, and you can select the desired care mode.

(3) Checking the Details of the Care and the Handpiece Effects


Check the details of the care and the handpiece effects, and then press the Next button.

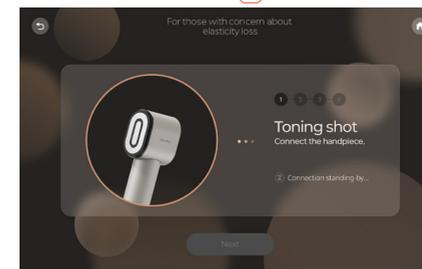(4) Connecting the Handpiece


The following handpiece connection window will appear, and the handpiece holder LED will turn on at the correct handpiece position.

2) How to Carry out Care in General Mode

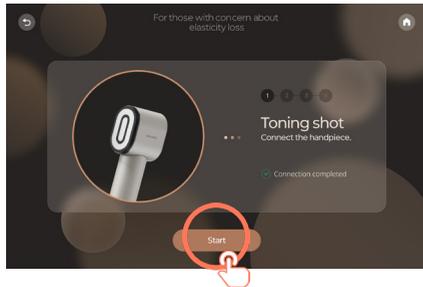(5) Selecting the Start Button


When the handpiece connection is completed, press the 'Start' button to start operation.

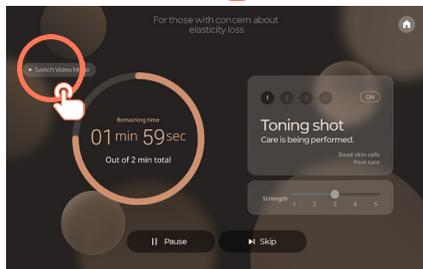(6) Carrying out the Selected Care Mode


The selected care will be carried out, and you can check the correct instructions for using the handpiece by clicking Video mode during operation.

※ For a detailed description of the screen, refer to page 14.

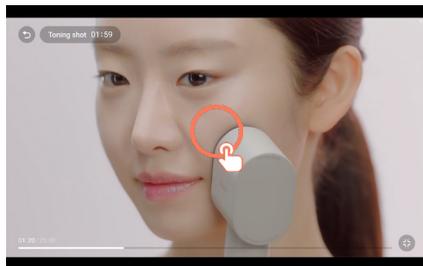(7) Options Settings


You can display the Pause, brightness settings and volume settings buttons by touching the video screen.

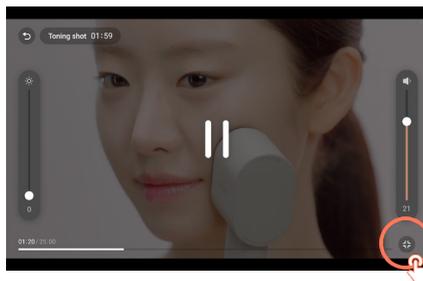(8) Pausing and Switching Screen Size

When you press the Pause button at the center, the video and the handpiece will be paused.


You can switch the mode to small screen mode by pressing the Minimize button at the bottom right side.

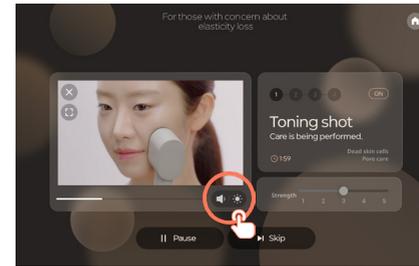(9) Small Screen Mode


In small screen mode, the timer will be displayed next to the clock icon on the right side of the video.

You can view the Settings pop-up window by pressing the Volume/Brightness Settings button.

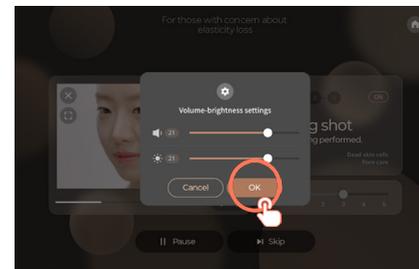(10) Volume/Brightness Settings


Adjust the volume and the brightness as you desire, and then press the OK button to return to the Care In Progress screen.

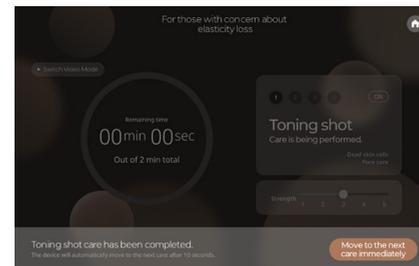(11) Carrying out Next Care Mode

When care is completed, the screen where you can press the 'Move to Next Care Immediately' button will appear.

If the button is not pressed, the screen will automatically move to the Next Care screen 10 seconds later.

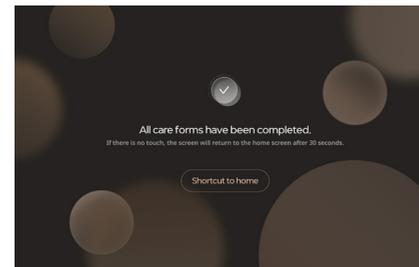
(12) Pause and Switching Screen Size

When care has completed for the last handpiece, the Completion screen will appear as follows.


If no button is pressed for 10 seconds, the screen will automatically move to the Home screen.

3) How to Use Expert Mode

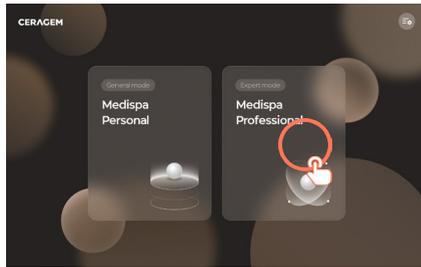
(1) Selecting Expert Mode

From the Home screen, select Expert mode.

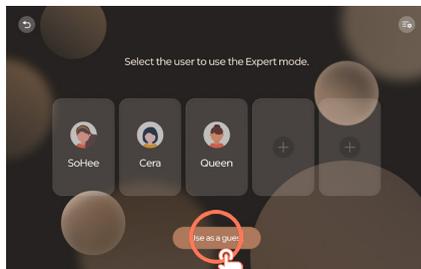(2) Selecting the User


If a user is not registered, you can click the 'Use as a guest' button to use Expert mode. When using as a guest, the Favorite function in Expert mode cannot be used.

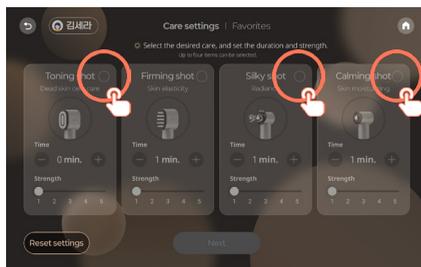(3) Selecting the Handpiece and Options

Select the desired handpiece and set the care time and operation strength.

※ For a detailed description of the screen, refer to page 15.

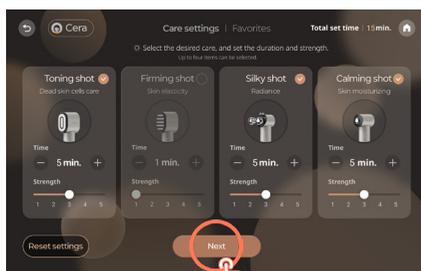
(4) Moving to the Screen

Press the Next button to move to the next screen.

(5) Starting after Checking the Details of Care


Check the details of the care you set, and then press Next to start care.

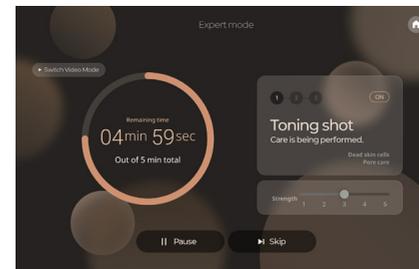(6) Carrying out Care

Care is carried out the same as in General mode.

In Expert mode, Video mode cannot be used.

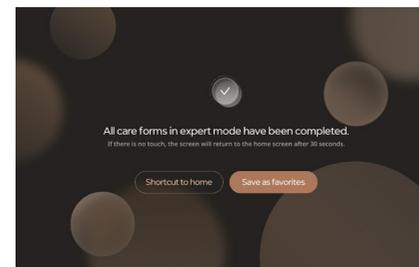
(7) Saving Favorites

When care has completed with the last handpiece, the Completion screen will appear as follows, and if no button is pressed for 30 seconds, the screen will automatically move to the Home screen.


In Expert mode, you can press the 'Save as Favorites' button to save the details of the care performed to 'my list of favorites' after the care has completed.

4) Favorites in Expert Mode

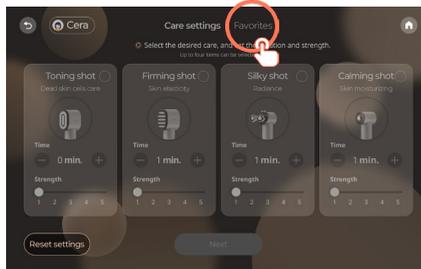
(1) Moving to the Favorites Tab

When logging in as a user, you can move to the “Favorites tab” in Expert mode and save the care with the desired details.

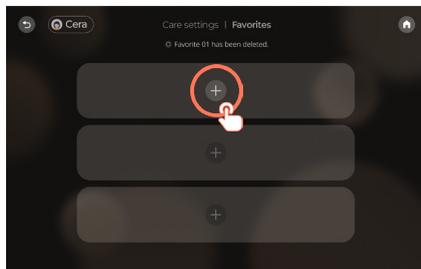
(2) Favorites Registration Settings Screen

When you press an empty space in the Favorites tab, the screen will move to the screen for setting the care to be registered.

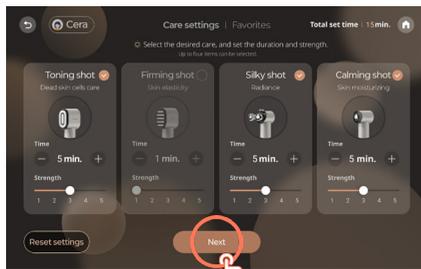(3) Saving Details of Favorites

When you set the details of the care to be registered and press the “Save’ button, the Favorites will be saved.

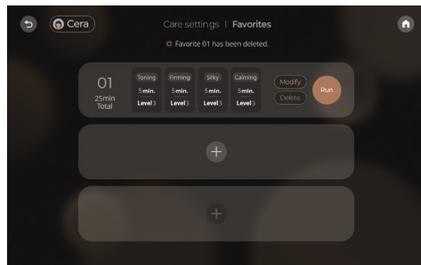
(4) List of Favorites

You can execute the favorite you registered from the Registered List of Favorites.

When you delete a favorite, you will not be able to recover it.


5) Configuration

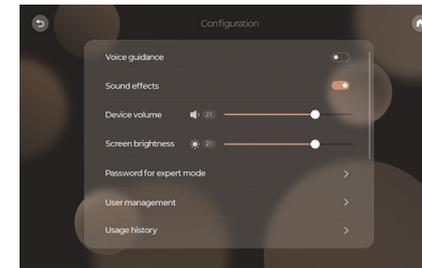
(1) Cellunic Configuration

You can check and change the following items from Configuration.

- Voice Guidance
- Sound Effects
- Device Volume
- Screen Brightness
- Password for Expert Mode
- User Management
- Usage History
- App Version, Etc.

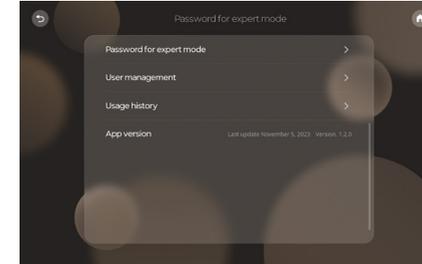
6) Configuration How to Save Usage History

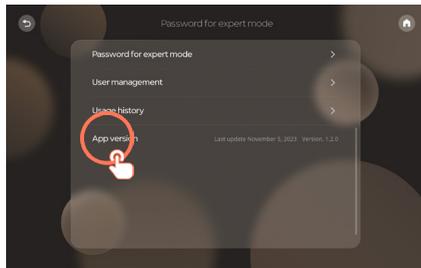(1) Checking the Usage History


From the Celluonic Configuration screen, press the Usage History item.

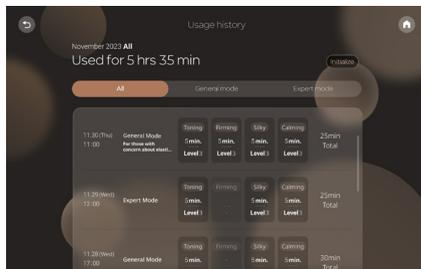(2) Usage History Display Screen


The 'Skipped' or 'Paused' time during operation is not displayed, and only the time for carrying out care and the strength of the handpiece at the end of usage is displayed.

※ Detailed information can be viewed when downloading Usage History.

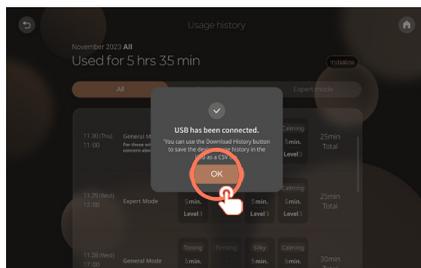(3) Connecting a USB device


When connecting a USB device to the USB port, the 'Download History' button will be created.

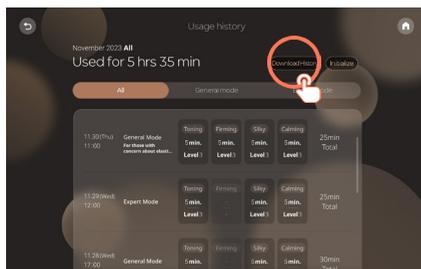(4) Downloading History


When you press the Download History button, the year, time and extension will be generated and saved automatically and these can be checked in the Excel file.

7) Configuration Initializing Usage History

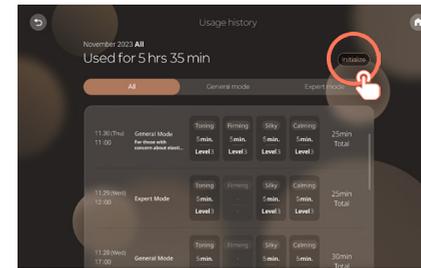(1) Usage History Screen

From the Usage History screen, press the Initialize button.


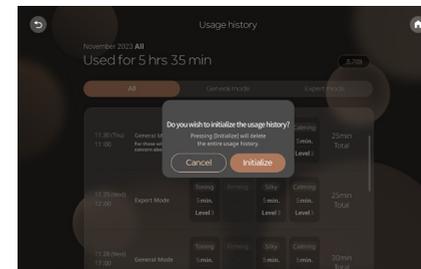
(2) Touching the Initialize Usage History Button

When you press the Initialize button on the displayed pop-up window, the Usage History will be initialized.

04. Troubleshooting

If there is a problem while using the product, first check the following before contacting the Customer Service Center (213-784-0475).

The product must not be disassembled or repaired arbitrarily by the user. If the product is damaged due to arbitrary disassembly or repair, the product warranty will not be recognized.

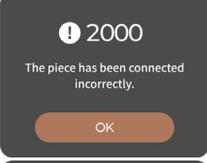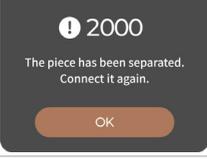Screen Display	Circumstances	Possible solutions
	When the handpiece cable is being wound	When operation is carried out after pressing the Wind Cable button, you cannot pull or wind it directly. Wait until the operation is completed.
	When the Wind button is pressed while the handpiece and the cable are connected	Detach the handpiece and the cable, and then press the Wind button.
	When the handpiece has been connected incorrectly while Care is set	Check the shape of the handpiece on the screen and connect the correct handpiece.
	When the handpiece has been detached while Care is being carried out	If the handpiece is detached while Care is being carried out, the operation will be suspended. If the handpiece has been separated, connect it again.

*If the symptom cannot be solved by referring to the suggestions herein, refer to the product warranty and contact the Customer Service Center.

05. FAQ

Q1. The device is operating, but I don't feel any stimulation.

Use the product after applying sufficient soothing gel or your preferred cosmetics. Detailed instructions for use can be checked on the Cellunic Website (ceragemus.com).

Q2. Can I use the product in the morning and evening?

Yes, you can. There may be individual differences, so please use the product according to your skin condition.

Q3. Can I use the same shot continuously?

Continuous use is not recommended, the product is set to be used for up to 10 minutes per single use.

Q4. Is it okay to use it every day?

It is recommended that you use it twice a week, and there may be individual differences, so please adjust the cycle of use, time and strength according to your skin condition.

Q5. How can I wash the product after use?

Wipe it off with a dry cotton pad or towel, and avoid using alcohol swabs or washing with water.

Q6. Can I select and use the desired mode?

You can use it in Expert mode. Use it according to the connection order of the handpiece.

06. Product Specification

Product Name	CERAGEM MEDISPA PRO		
Model Name	CGM ESK-2405		
Product Dimension	11.3 (W) X 8.3 (L) X 16.0 (H) in		
FEPS	CLASS III EQUIPMENT (Form of electric Shock protection)		
Product Weight	Entire Device	9.3 lbs	
	Main Body	7.9 lbs	
	Handpiece	Toning Shot 0.33 lbs	Firming Shot 0.30 lbs
		Silky Shot 0.45 lbs	Calming Shot 0.23 lbs
Input Power	100-240V~, 50/60Hz		
Output Power	DC 12V ---, 5A		
NRTL No.	E116039 Complies with UL 1431 CSA C22.2 No. 60335-1		
FCC No.	KR24-YEF0110		
CS Center	213-784-0475 (Customer Service)		
Distributor	Ceragem International Inc. 9350 Wilshire Blvd. Suite 400 Beverly Hills, CA 90212 info@ceragemus.com		
Manufacturer	Easytem Co., Ltd. 25 Podowon-ro 116 beon-gil, room 304, Siheung-si, Gyeonggi-do, Republic of Korea		

* This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions :
(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

07. Product Warranty

1. Ceragem provides warranty services for this product in accordance with the regulations under Consumer Protection Law.
2. The warranty period is calculated from the date of purchase, so have the exact date of purchase written on the product warranty.
3. If you request a repair due to product failure, contact us or the Customer Service Center at 213-784-0475.

Product Warranty	Product Name	Model Name
	Serial No.	Place of Purchase
	Date of Purchase	Date of installation
Customer Information	Customer's Name	
	Telephone Number	
	Address	
	Location	
	Zip Code	
	※ Warranty Period: It is 1 year (12 months) from the date of purchase, or in the case of rental, it varies according to the agreed rental period.	
	※ Parts Holding Period: Consumer dispute resolution standards apply.	

- Information**
- (1) The warranty period is 1 year from the date of purchase, and if it is difficult to confirm the date of purchase, 15 months from the date of manufacture apply as the warranty period.
 - (2) During the warranty period, any product failures or damages related to the following circumstances are not covered by this warranty.
 - Maintenance and repair performed by persons not approved by the manufacturer
 - When damage occurred to the product due to another company's product other than the product provided by the manufacturer
 - When damage occurred to the product due to the use of a part or a consumable not specified by the manufacturer
 - Use of an accessory not specified by the manufacturer
 - Damage due to an accident or a natural disaster (earthquake, flooding, etc.)
 - Negligence by the consumer
(product failure due to misuse of the product or incorrect repair by the consumer)
 - Consumable parts reaching the end of their life span (battery, etc.)
 - Damage and normal wear, minor defects that do not interfere with use
 - Damage caused by failure to comply with the precautions and operating instructions in the user manual
 - Damage caused by ignoring specified maintenance and inspection

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